

M



- World class German engineering
- Cobalt Chrome, LINK PorEx¹ and cementless options
- Wide range of sizes 14 femoral & 10 tibial
- 2up 2down femoral tibial size compatibility²
- CR, PS, PS+, All-Poly bearings in UHMWPE & Vitamin E infused cross-linked polyethylene with 1mm increments³
- Elegant, reproducible instruments configured in 3 trays
- Cost conscious implant choices ideal for same-day surgery

¹Titanium Niobium Nitride Coating. ²2 up/down not available in all sizes. ³Imm increments for certain implant constructs only. ⁴Internal Technical Report.



- Complete revision system suitable for standard and complex anatomical situations
- FlexiCones friendly tibial design
- Vast array of straight and offset cemented and cementless stems
- Precise instrumentation streamlined into four trays¹

¹Most cases performed with 4 trays.


Soft Tissue Friendly Anterior Flange


> Wide Proximal Trochlea

Harmonic Growth Between Sizes (3mm)

> Lateralized Trochlear Groove

Anatomical Notch Geometry

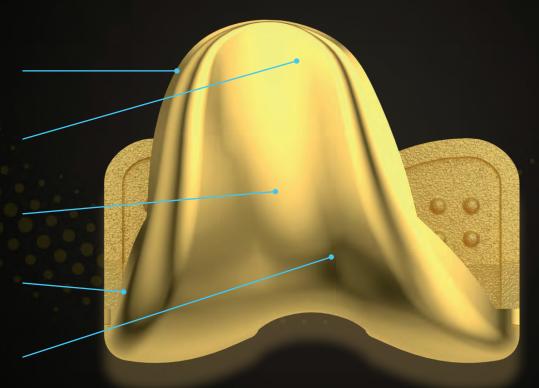

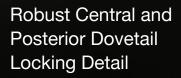

Antirotation Island Guides Bearing Implantation

Flattened Anterior Lip Designed To Reduce Bearing Insertion Angle

Keel and Fin Designed to Protect Against Rotational Forces 4.5° Anterior Flange

Low Profile Patella Flange

Gently Decreasing Flexion Radius

> Graduated Trochlear Groove Depth

> > LSK CR Femur

Bone Conserving Box Geometry

Post Engagement @ 60° Flexion

High PS Cam Location

Spherogrip Macrostructure Provides for a Homogeneous Bone Cement-Implant Interface And Uniform Force Transmission in the Cement Mantle⁴

LSK PS Femur



Manufactured in Germany



This advertising is not intended to minimize the risks associated with surgery and implantation, nor does it contain any claim as to individual risk profiles of patients, operation schedules, or therapeutic decisions of surgeons. In particular, but not limited to, it does not mean to override any product related safety requirements. Instructions for use for any product manufactured and distributed by the Waldemar Link GmbH & Co. KG, Germany, apply without limitation. In case of any doubt, please contact our US subsidiary LinkBio for further information.

Waldemar Link GmbH & Co. KG and/or other corporate affiliated entities own, use or have applied for the following trademarks in many jurisdictions: LINK, BiMobile, SPII, Modell Lubinus, E-Dur, EndoDur, T.O.P. II, BetaCup, CombiCup PF, CombiCup SC, CombiCup R, MobileLink, C.F.P., LCU, SP-CL, LCP, MIT-H, Endo-Model, Endo-Model SL, MP, MEGASYSTEM-C, GEMINI SL, SPAR-K, LCK, Link OptiStem, HX, TiCaP, X-LINKed, PorAg, LINK PorEx, BiPorEx, PorEx-Z, TrabecuLink, Tilastan, customLINK, RescueSleeve, Stactip, VACUCAST.

Other trademarks and trade names may be used in this document to refer to either the entities claiming the marks and/or names or their products and are the property of their respective owners.

Waldemar Link GmbH & Co. KG

Barkhausenweg 10 · 22339 Hamburg · Germany Phone +49 40 53995-0 · info@link-ortho.com www.link-ortho.com

